


**ORPORT ICU** pendant is designed for medical equipment integration. Base on clinical needs, ORPORT ICU pendant combines diagnostic equipment, patient monitor, medical ventilator, IV pole, injection pump, gas outlets and power sockets.

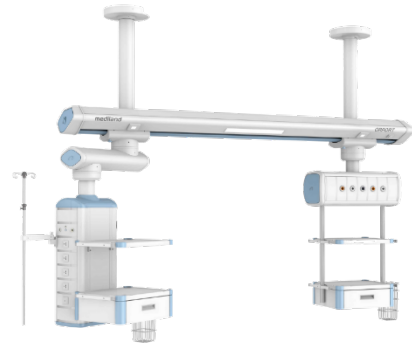

### Key Features of ORPORT ICU:

#### Stability and Flexibility

Built from single-molded aluminum alloy, it is able to withstand heavy equipment. The linear slide rail design also allows smooth and easy horizontal movements.

#### Soft LED lighting effect

Multiple LED light sources on the crossbeam produce a soft light environment for medical staff.

#### Extendable port

Consistent port design gives customers the choice to extend various gas outlets, power socket and communication/ image ports for future demand.

## Accessories

Drawer, shelf, IV pole, stainless steel basket, monitor rack, examination light, ambient light and cable clamp can be purchased optionally.

## Customized design

A variety of combination of service heads and arms are available based on customer needs.

## Technical Data

Material	Single-molded aluminum alloy
Weight Capacity of Crossbeam	Max.350kg
Weight Capacity of Single Pendant	Max.150kg
Length of Crossbeam	$\geq 2100$ mm
Horizontal moving Diameter	500mm
Brake System	Standard damper braking system